

Performance of an in-vivo tear film osmometer in normal ocular surface conditions.

Dr. Henry Reis, MD; Stefanie Grenier, BSc; Dr. Daniela Albuquerque, MD.

Purpose

To assess the performance of an in-vivo osmometer in a clinical setting using the i-Pen osmometer in patients who scored within normal limits on the SPEED questionnaire.

Methods

A total of 28 patients (ages 26-87) underwent comprehensive ocular surface assessment at Focus Eyecare Centre, an accredited dry eye centre in Burnaby, BC, Canada. Patients were screened for ocular surface disease (OSD) with the SPEED questionnaire and the in-vivo osmometer i-Pen, and included in the study if all eligibility criteria were met.

Ocular Surface Disease (OSD) has been objectively correlated to elevated tear osmolarity, and most progressive diagnostic algorithms currently include osmometers. Eye care practitioners must rely on objective testing to better evaluate patients with dry eye disease, whose gamut of symptoms covers a wide range of presentations from mild ocular pruritus to severe ocular pain.

A fundamental concern in OSD investigation is that the performance of osmometers could be influenced by a variety of factors. Among those that could lead to osmolarity underestimation are the administration of non-contact tonometry, retinal photography and topical medication up to two hours prior to OSD screening. Other factors lead to osmolarity overestimation, such as the use of eye makeup and/or contact lenses, or even the excessive use of handheld electronic devices on the day of the examination.

It is paramount to understand that a test with high specificity allows eye care practitioners to rule out a condition when encountering a negative result. Several studies have been conducted to assess the validity of osmolarity testing, but most study designs include patients with mild, moderate or severe disease, with emphasis in test sensitivity. Contrary to that approach, the authors intended to evaluate the true negative rate (TNR), also known as test specificity (SPC) of the commercially available in-vivo osmometer, i-Med Pharma's i-Pen, in patients with normal ocular surface.

$$TNR = \frac{\sum TN}{\sum \text{condition negative}}$$

The study included 56 eyes, initially screened with the SPEED questionnaire. Patients who failed to meet eligibility criteria were excluded from the study. Inclusion criteria included: patient age greater than 18 years, negative history of makeup or topical ocular medication within 24 hours from the exam, negative history of ocular disease, and SPEED score equal or inferior to 8. Patients were then subcategorized by SPEED score, age, ethnicity and time of osmolarity measurement.

Osmolarity measurements with the i-Pen were acquired by CCOA (Canadian Certified Optometric Assistant) personnel strictly following the manufacturer's manual, with proper patient consent.

Results

According to the i-Pen user manual, the mean osmolarity value in normal patients was 300 mOsm/L and a range from 275 to 316 mOsm/L was documented.

The mean tear osmolarity encountered in this study was 295 mOsm/L, with an average of 300.85 mOsm/L. Exam values ranged from 276 to 336 mOsm/L, with a statistically higher upper reference limit (URL) when compared to the URL published by the manufacturer. In view of that discrepancy, one might claim that most symptoms of OSD have no defined, sudden onset, leading to a proportionally small group of those patients being classified as normal. OSD symptoms might be initially masked or subclinical, leading to artificially low questionnaire scores in patients that would otherwise be classified as positive for ocular surface disease.


The analysis of tear osmolarity in different ethnical groups unveiled a statistically significant variance, with patients of African American descent presenting with the lowest average value (287 mOsm/L), followed by Mixed (292 mOsm/L), Caucasians (298.36 mOsm/L), Hispanic (302.22 mOsm/L) and lastly Asians (312.8 mOsm/L).

No statistical significance was observed when analyzing the distribution of tear osmolarity by age or time of testing.


The i-Pen True Negative Rate (TNR), also known as Specificity (SPC), was calculated as 91.7%. This is considered clinically significant, as most tests that present SPC values over 90% are considered reliable; despite the fact that a quintessential gold standard test would display both sensitivity and specificity of 100%. Furthermore, a direct TNR comparison between the two commercially available osmometers, Tearlab and i-Pen indicates that both options have similar test specificity, with 88% for the former and 91.7% for the latter.

Conclusions


The inclusion of tear osmolarity in ocular surface disease diagnostic algorithms is strongly encouraged by the authors. The test has proved to have high specificity, allowing eye care practitioners to exclude OSD in most cases where osmolarity testing is negative. It is also important to consider the statistically significant correlation between ethnicity and ocular surface osmolarity, and be conscientious when analyzing dry eye questionnaires, as some patients might neglect to report symptoms in early stages of the condition.

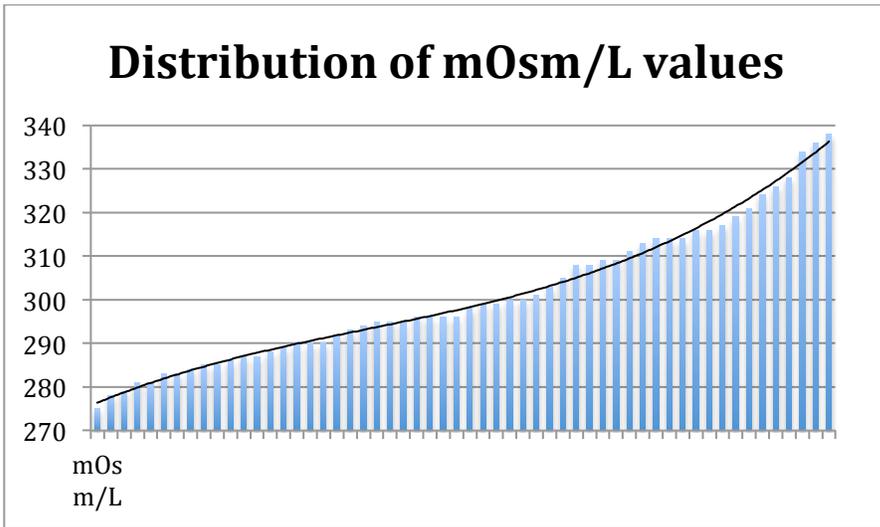

Distribution of tear osmolarity in patients with normal ocular surface conditions.

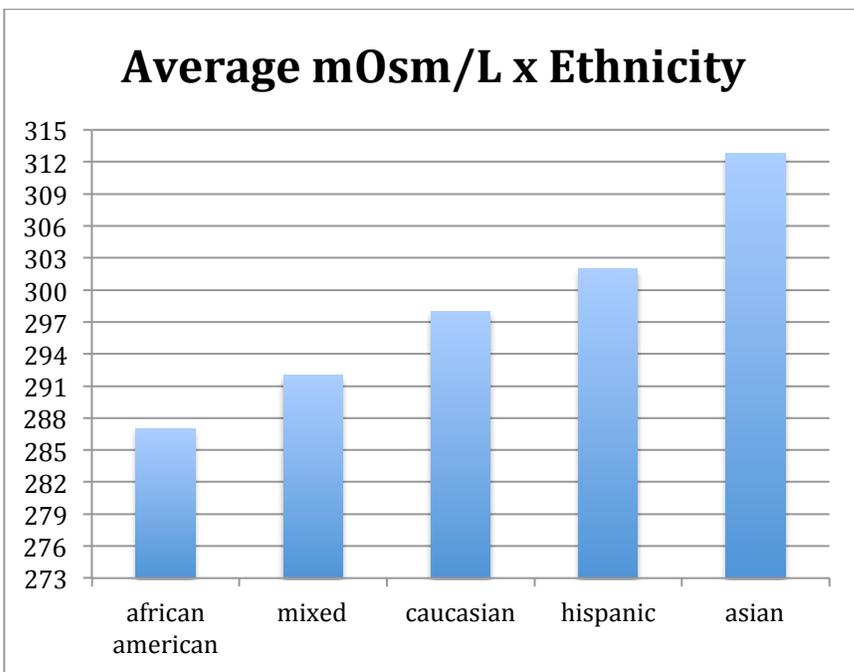

Analysis of tear osmolarity versus ethnicity

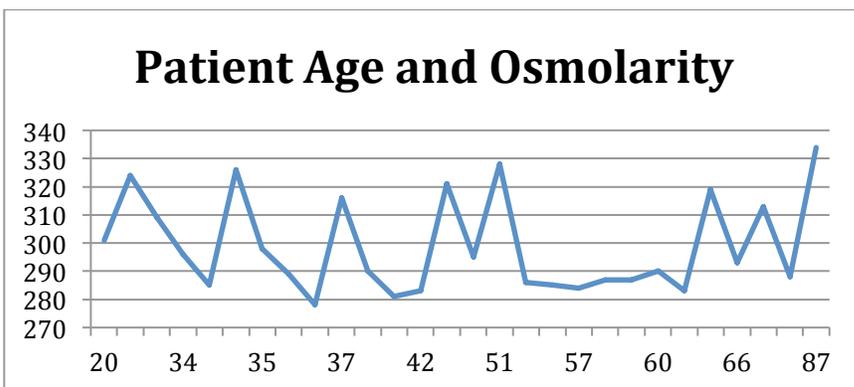

Analysis of tear osmolarity versus patient age yielded no statistical significance.

1. Karpecki PM. Diagnosis and treatment of ocular surface conditions: focus on dry eyes. *Clinical and Refractive Optometry* 2016; 27(6): 236-241.
2. Definition and Classification of Dry Eye. Report of the Diagnosis and Classification Subcommittee of the Dry Eye Workshop (DEWS). *The Ocular Surface* 5(2): 75-92, 2007.
3. Lemp MA, Foulks GN. The Definition & Classification of Dry Eye Disease Guidelines from the 2007 International Dry Eye Workshop. *Ophthalmology Management* April 2008.
4. Tomlinson A, Khanal S, Ramaesh K, et al. Tear film osmolarity: determination of a referent for dry eye diagnosis. *Invest Ophthalmol Vis Sci.* 2006 Oct;47(10):4309-15.
5. Rocha G, Gulliver R, Borovik A, Chan CC. Randomized, masked, in vitro comparison of three commercially available tear film osmometers. *Clinical Ophthalmology* 2017; 11: 243-248.
6. Versura P, Profazio V, Campos EC. Performance of tear osmolarity compared to previous diagnostic tests for dry eye diseases. *Current Eye Research* 2010; 35(7): 553-564.
7. Keir N, Ngo W, Situ P, Korb D, et al. Evaluation of the standard patient evaluation of eye dryness (SPEED) questionnaire. *Investigative Ophthalmology & Visual Sciences.* 2013; 54: 6028.
8. Ngo W, Srinivasan S, Keech A, Keir N, et al. Self versus examiner administration of the ocular surface disease index. *J Optom.* 2017; 10(1): 34-42.
9. Narayanan N. Osmolarity: a diagnostic test for dry eyes. *Review of Optometry* 2011.
10. Tauber S. Treating and managing dry eye. *Clinical and Refractive Optometry* 2016; 27(1): 14-17.
11. Avinashi P. Principles of dry eye disease: diagnosis, treatment and management. *Clinical and Refractive Optometry* 2016; 27(4): 127-132.
12. Barrus BJ. The science of dry eye hyperosmolarity. *Clinical and Refractive Optometry* 2016; 27(5): 202-206.
13. Maharaj RL. In vivo ocular surface osmolarity in a dry eye population. *Clinical and Refractive Optometry* 2017; 28(1): 3-6.
14. Baudouin C, et al. Role of hyper- osmolarity in the pathogenesis and management of dry eye disease. *Ocul Surf* 2013; 11: 246-258.
15. Keech A, Jones L, Senchyna M. TearLab Osmometer Validation Protocol. Centre for Contact Lens Research. University of Waterloo School of Optometry. March 2009.
16. I-Pen Osmolarity System User Manual. I-Med Pharma Inc, 2016.